

K971448

## **SUMMARY OF SAFETY EFFECTIVENESS**

**000042**

JUL 3 1997

### **I. General Provisions**

Common or Usual Name: PTA Balloon Catheter

Proprietary Name: Opta LP™ PTA Balloon Catheter

### **II. Name of Predicate Devices**

- Cordis Opta5 PTA Balloon Catheter
- Cordis PowerFlex PTA Balloon Catheter
- Cordis Small Vessel PTA Balloon Catheter
- Cordis Savvy PTA Balloon Catheter
- Meditech UltraThin Balloon Catheter

### **III. Classification**

Class II

### **IV. Performance Standards**

Performance standards have not been established by the FDA under section 514 of the Food, Drug and Cosmetic Act.

### **V. Intended Use and Device Description**

The Opta LP PTA balloon catheters is intended to dilate stenoses in iliac, femoral, ilio-femoral, popliteal infra popliteal and renal arteries, and for treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae.

The Cordis Opta LP PTA balloon catheter is a dual lumen design with a distal inflatable balloon. Two radiopaque marker bands indicate the dilating section of the balloon and aid in balloon placement.

The working pressure range for the balloon is between the nominal size pressure and the rated burst pressure. All balloons distend to sizes above the nominal size at pressures greater than the nominal pressure.

The balloon inflation lumen is used to inflate and deflate the balloon. The nominal balloon size is printed on the hub.

The guidewire lumen is used to track the catheter over a prepositioned guidewire or to inject contrast medium and/or saline. The maximum injection pressure is 450 psi. The compatible guidewire size, catheter shaft French size and catheter length are printed on the hub. The radiopaque marker bands indicate the stated nominal length of the balloon.



**VI. Biocompatibility**

All materials used in the Opta LP PTA balloon catheter are biocompatible.

**VII. Summary of Substantial Equivalence**

The Cordis Opta LP PTA balloon catheter and the referenced Cordis Opta5 PTA balloon catheter, Cordis PowerFlex PTA balloon catheter, Cordis Small Vessel PTA balloon catheter, Cordis Savvy PTA balloon catheter and Meditech UltraThin balloon catheter are similar in their basic design, construction, indication for use and performance characteristics.





DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

Mirjam Barboza, M.D.  
Manager, Regulatory Affairs  
and Clinical Research  
Cordis Corporation  
P.O. Box 025700  
Miami, Florida 33102-5700

JUL 3 1997

Re: K971448  
Cordis Opta™ LP PTA Balloon Catheter  
Regulatory Class: II (two)  
Product Code: LIT  
Dated: April 18, 1997  
Received: April 21, 1997

Dear Dr. Barboza:

We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.


If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the Current Good Manufacturing Practice requirement, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic QS inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.



This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices); please contact the Office of Compliance at (301) 594-4648. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its internet address "<http://www.fda.gov/cdrh/dsmamain.html>".

Sincerely yours,

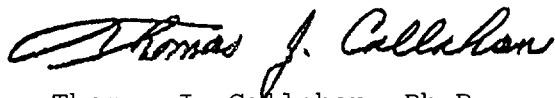A handwritten signature in black ink, reading "Thomas J. Callahan". The signature is fluid and cursive, with the first name "Thomas" being the most prominent part.

Thomas J. Callahan, Ph.D.  
Director  
Division of Cardiovascular, Respiratory,  
and Neurological Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure




000060

510(k) Number (if known):

Device Name: Cordis Opta LP™ PTA Balloon Catheter

Indications for Use:

**The Opta LP PTA Balloon Catheter is intended to dilate stenoses in iliac, femoral, ilio-femoral, popliteal, infra popliteal and renal arteries, and for treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae.**

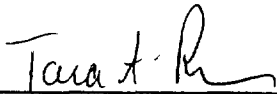

(Division Sign-Off)

Division of Cardiovascular, Respiratory,  
and Neurological Devices

510(k) Number K971448

PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IS NEEDED

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ✓  
(Per 21 CFR 801.109)

OR

Over-The-Counter Use \_\_\_\_\_