

5/5/99

## Attachment D

# 510(k) Summary

The assigned 510(k) number is K991249

### Submitter Information (21 CFR 807.92(a)(1))

Submitter: Becton Dickinson  
1 Becton Drive  
Franklin Lakes, NJ 07417-1880

Contact: Cindy Morrow  
Sr. Regulatory Specialist  
(408) 954-2694  
(408) 954-2495 (FAX)  
cmorrow@bdis.com

Summary date: April 12, 1999

### Device Name/Classification (21 CFR 807.92(a)(2))

Name: FALCON® IVF Four Well Plate

Classification: Assisted reproduction labware, Class II, 884.6160, Code: 85 MQK

### Substantially Equivalent/Predicate Device (21 CFR 807.92(a)(3))

This product is being submitted according to the Federal Register notice located on page 48428 of Vol.63, No. 175 on September 10, 1998 under the title Obstetric and Gynecologic Devices; Reclassification and Classification of Medical Devices Used for In Vitro Fertilization and Related Assisted Reproduction Procedures.

### Device Description (21 CFR 807.92(a)(4))

The FALCON® IVF Four Well Plate is sterile (SAL of  $10^{-6}$ ), non-pyrogenic by Limulus Amebocyte Lysate (LAL of  $< 20$  EU/device), and nonembryotoxic as tested by the mouse embryotoxicity assay (MEA) 2-cell method. The single-use plastic plate has four wells, each well area is 1.39 cm<sup>2</sup> and each well volume is 1.8 mL. The plate has a unique lid which provides access to two wells at a time, while two remain covered. The wells are numbered and a large writing patch allows clear sample identification. Plates are packaged in individual peel-open trays for sterile presentation and are shipped in cases of 100 plates per case.

The plates have perfectly flat, optically clear surfaces for optimum manipulation and observation of the ova and embryos. The lids are designed for aseptic manipulation and


consistent venting to maintain proper humidification. The plates are manufactured from virgin crystalline polystyrene tested for USP Class IV, V, and VI cytotoxicity. and the surfaces are treated to provide a more wettable or hydrophilic surface for tissue culture.

**Intended Use (21 CFR 807.92(a)(5))**

The FALCON® IVF Four Well Plate is intended for use in preparing, storing, manipulating, or transferring human gametes or embryos for in vitro fertilization (IVF), gamete intrafallopian transfer (GIFT), or other reproduction procedures.



David Ball  
Director of Quality Assurance/ Regulatory Affairs  
Becton Dickinson

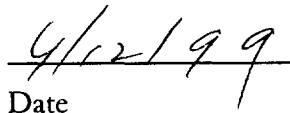

Date



Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

MAY 5 1999

Mr. David Ball  
Director of Quality Assurance/ Regulatory Affairs  
Becton Dickinson & Company  
1 Becton Drive  
Franklin Lakes, NJ 07417-1884

Re: K991249  
FALCON® IVF Four Well Plate  
Dated: April 12, 1999  
Received: April 13, 1999  
Regulatory Class: II  
21 CFR 884.6160/Procode: 85 MQK

Dear Mr. Ball:

We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the Current Good Manufacturing Practice requirements, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic QS inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4613. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597, or at its internet address "<http://www.fda.gov/cdrh/dsma/dsmamain.html>".

Sincerely yours,

CAPT Daniel G. Schultz, M.D.  
Acting Director, Division of Reproductive,  
Abdominal, Ear, Nose and Throat,  
and Radiological Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

**Attachment C**

# Indications for Use

510(k) Number: K 991249

Device Name: FALCON<sup>®</sup> IVF Four Well Plate

The FALCON<sup>®</sup> IVF Four Well Plate is sterile, nonpyrogenic, embryotoxicity tested, single-use plasticware intended to prepare, store, manipulate, or transfer human gametes or embryos for in vitro fertilization (IVF) or other assisted reproduction techniques.

(Please do not write below this line – Continue on another page if needed)

---

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ✓  
(Per 21 CFR § 801.109)

Or

Over-the-Counter Use ~~XXXXXX~~

510(k) Premarket Notification

  
(Division Sign-Off)

Division of Reproductive, Abdominal, ENT,  
and Radiological Devices

510(k) Number

K991249