

DEC 4 1998

K983648
EXHIBIT #1
Page 1 of 2

510 (K) SUMMARY

This summary of 510(k) safety and effectiveness information is being submitted in accordance with the requirements of SMDA 1990 and 21 CFR §807.92.

The assigned 510(k) number is: 983645

1. Submitter's Identification:

Mr. Fang Yi Liu
Shijiazhuang Great Eagle
Plastic Products Co., Ltd.
685 He Ping Road
New Development Zone
Shijiazhuang, He Bei Province
P.R. China

Date Summary Prepared: October 14, 1998

2. Name of the Device:

Shijiazhuang Great Eagle Plastic Products Co., Ltd.
Powdered Vinyl Patient Examination Gloves

3. Predicate Device Information:

Sunmax Enterprise Shanghai Co., Ltd. Powdered Vinyl Patient
Examination Gloves, K#960746
Shanghai Super Gloves Company, Ltd. Powdered Vinyl Patient
Examination Gloves, K#974151

4. Device Description:

Classified by FDA's General and Plastic Surgery Device Panel
as Class I, 21 CFR 880.6250, Powdered Vinyl Patient
Examination Glove, 80LYZ, Powdered with Absorbable Dusting
Powder, USP, Class III and meets all requirements of ASTM
Standard D5250-92.

5. Intended Use:

A patient examination glove is a disposable device intended
for medical purposes that is worn on the hand of healthcare
and similar personnel to prevent contamination between
healthcare personnel and the patient's body, fluids, waste or
environment.

6. Comparison to Predicate Devices:

Shijiazhuang Great Eagle Plastic Products Co., Ltd. Synthetic Powdered Vinyl Patient Examination Gloves, is substantially equivalent in safety and effectiveness to the Sunmax Enterprise Shanghai Co., Ltd. Powdered Vinyl Patient Examination Gloves and the Shanghai Super Gloves Company, Ltd. Powdered Vinyl Patient Examination Gloves.

7. Discussion of Non-Clinical Tests Performed for Determination of Substantial Equivalence are as follows:

The standards used for Shijiazhuang Great Eagle Plastic Products Co., Ltd. glove production are based on ASTM-D-5250-92. All testing meets requirements for Physical and Dimensions Testing conducted on gloves, Inspection Level S-2, AQL 4.0.

The FDA 1000 ml. Water Fill Test was also conducted with samplings of AQL 2.5, Inspection Level S-4, meeting these requirements. Primary Skin Irritation and Skin Sensitization (allergic contact dermatitis) testing was conducted with results showing no primary skin irritant or sensitization reactions.

There are no special labeling claims and we do not claim our gloves as hypoallergenic on our labels.

8. Discussion of Clinical Tests Performed:

Not Applicable - There is no hypoallergenic claim.

9. Conclusions:

Shijiazhuang Great Eagle Plastic Products Co., Ltd. Synthetic Vinyl Patient Examination Gloves (Powdered) conform fully to ASTM-D-5250-92 standards as well as applicable 21 CFR references, and, meets pinhole FDA requirements, biocompatibility requirements and labeling claims as shown by data in Section 7. There are no safety/efficacy issues or new claims from the "substantial equivalence" products cited.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

DEC 4 1998

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

Shijiazhuang Great Eagle Plastics Products
Company, Ltd.
C/O Ms. Susan Goldstein-Falk
Official Correspondent
Basic Medical Industries, Incorporated
13945 Magnolia Avenue
Chino, California 91710

Re: K983645
Trade Name: Powdered Vinyl Patient Examination Gloves
Regulatory Class: I
Product Code: LYZ
Dated: October 14, 1998
Received: October 16, 1998

Dear Ms. Goldstein-Falk:

We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.


If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the Good Manufacturing Practice for Medical Devices: General (GMP) regulation (21 CFR Part 820) and that, through periodic GMP inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

Page 2 - Ms. Goldstein-Falk

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4692. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its internet address "<http://www.fda.gov/cdrh/dsmamain.html>".

Sincerely yours,

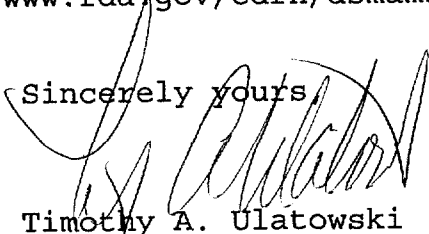

Timothy A. Ulatowski
Director
Division of Dental, Infection Control,
and General Hospital Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

510(k) NUMBER (IF KNOWN): _____

DEVICE NAME: Shijiazhuang Great Eagle Plastics Products Co., Ltd.INDICATIONS FOR USE: Powdered Vinyl Patient Examination Gloves-Powdered

A patient examination glove is a disposable device intended for medical purposes that is worn on the examiner's hand or finger to prevent contamination between patient and examiner.

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE
IF NEEDED.)

Concurrent of CDRH, Office of Device Evaluation (ODE)

Prescription Use
(Per 21 CFR 801.109)

Chin S. Lin
(Division Sign-Off)

Division of Dental, Infection Control,
and General Hospital Devices

OR Over-The-Counter-Use X
(Optional Format 1-2-96)

510(k) Number K 983645