

STERIS®



**510(k) Summary
For
SYSTEM 1E Liquid Chemical Sterilant Processing System**

STERIS Corporation
5960 Heisley Road
Mentor, OH 44060
Phone: (440) 354-2600
Fax No: (440) 357-9198

Contact: Robert Sullivan
Senior Director, FDA Regulatory Affairs
Tel: 440-392-7695
Fax: 440-357-9198

Submission Date: November 28, 2011

1. Device Name

Trade Name: **SYSTEM 1E Liquid Chemical Sterilant Processing System**

Common/usual Name: Liquid Chemical Sterilizer

Classification Name: Sterilant, Medical devices, Liquid Chemical Sterilants/Disinfectants

Classification Number: 21 CFR 880.6885

Product Code: MED

2. Predicate Device

SYSTEM 1E Liquid Chemical Sterilant Processing System
K090036 modified by K101409, K102462 and addition to file

3. Description of Device

The SYSTEM 1E Liquid Chemical Sterilant Processing System is a liquid chemical sterilization system, utilizing peracetic acid to process totally immersible, heat sensitive, flexible and rigid endoscopes and their accessories, and microsurgical instruments. The system consists of the SYSTEM 1E Processor and the S40 Sterilant Concentrate, interchangeable processing trays/containers and Quick Connects. The device was originally cleared under K090036. Following this initial clearance the FDA was notified of a second qualified supplier of an ultraviolet light water treatment subsystem providing equivalent water treatment performance, safety, and effectiveness to that described within K090036. This modification was subsequently cleared under K101409. Software updates were cleared in K102462 and additions to that file.

The SYSTEM 1E Processor is an automated, self-contained device which creates and maintains the conditions necessary for liquid chemical sterilization in 6 minutes. Following the processing, the liquid chemically sterilized articles are rinsed with extensively treated water produced by passing EPA potable tap water through pre-filters, an ultraviolet light treatment subsystem, and then through two 0.1 micron filter membranes. The processor, which is computer controlled and continually monitored, provides printed documentation of each cycle.

S40 Sterilant Concentrate is a single use chemical sterilant concentrate developed exclusively for use in the SYSTEM 1E Processor. The active ingredient in S40, peracetic acid, is combined with inert ingredients to form a use dilution which inhibits corrosion of metals, polymers and other materials.

The interchangeable processing trays/containers are made to accommodate a variety of instrument types, models and procedure specific sets. Each container is designed to maintain instruments in appropriate position while specific Quick Connects ensure thorough circulation of the sterilant solution and rinses.

4. **Intended Use**

The SYSTEM 1E Liquid Chemical Sterilant Processing System is intended for liquid chemical sterilization of cleaned, immersible, and reusable critical and semi-critical heat-sensitive medical devices in healthcare facilities.

The SYSTEM 1E Processor dilutes the S40 Sterilant Concentrate to its use dilution (>1820 mg/L peracetic acid), liquid chemically sterilizes the load during a controlled 6-minute exposure at 46 to 55°C, and rinses the load with extensively treated* potable water. The processed load should be used immediately.

The SYSTEM 1E Processor uses only S40 Sterilant Concentrate to liquid chemically sterilize medical devices.

- * The extensive treatment of EPA potable water consists of:
 1. Pre-filtration through two pre-filters:
 - Pre-filter 1 is a gross depth filter that removes approximately 5 micron or larger particles/contaminants.
 - Pre-filter 2 is a surface filter that removes particles/contaminants > 0.1 micron.
 2. UV Irradiation:
 - During transit through the UV water treatment chamber, a UV dose sufficient to achieve a ≥ 6 -log reduction of MS2 virus is delivered to the water.
 3. 0.1 micron filtration:
 - The water prepared by pre-filtration and UV irradiation is filtered through redundant, 0.1-micron (absolute rated) membranes to remove bacteria, fungi and protozoa > 0.1 micron.

5. **Description of Safety and Substantial Equivalence**

The SYSTEM 1E Liquid Chemical Sterilant Processing System is the same as the predicates (K090036, K101409 and K102462) except for the modifications to specific components and software described in this submission and corresponding necessary changes to the Operating Manual for consistency. These modifications are intended to improve device usability, reliability, and serviceability without affecting any critical performance parameters. The modifications raise no new concerns of safety and effectiveness when compared to the predicate device.

The following table summarizes the verification activities that were performed, with their respective acceptance criteria, to ensure that the software modifications do not affect the safety or effectiveness of the SYSTEM 1E Liquid Chemical Sterilant Processing System.

Test	Acceptance Criteria	Result
Software Verification /Code Walk-Through	Only targeted modifications were implemented to the device software.	Pass
Software Validation	Only targeted modifications were implemented to the device software and modifications work as intended.	Pass



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
10903 New Hampshire Avenue
Document Control Room - WO66-G609
Silver Spring, MD 20993-0002

Mr. Robert Sullivan
Senior Director, FDA Regulatory Affairs
STERIS Corporation
5960 Heisley Road
Mentor, Ohio 44060

FEB 21 2012

Re: K113520
Trade/Device Name: SYSTEM 1E Liquid Chemical Sterilant Processing System
Regulation Number: 21 CFR 880.6885
Regulation Name: Liquid Chemical Sterilants / High Level Disinfectants
Regulatory Class: II
Product Code: MED
Dated: January 27, 2012
Received: January 30, 2012

Dear Mr. Sullivan:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration. Please note: CDRH does not evaluate information related to contract liability warranties. We remind you, however, that device labeling must be truthful and not misleading.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); medical device reporting (reporting of medical device-related adverse events) (21 CFR 803); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please go to <http://www.fda.gov/AboutFDA/CentersOffices/CDRH/CDRHOffices/ucml15809.htm> for the Center for Devices and Radiological Health's (CDRH's) Office of Compliance. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). For questions regarding the reporting of adverse events under the MDR regulation (21 CFR Part 803), please go to <http://www.fda.gov/MedicalDevices/Safety/ReportaProblem/default.htm> for the CDRH's Office of Surveillance and Biometrics/Division of Postmarket Surveillance.

You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 796-7100 or at its Internet address <http://www.fda.gov/MedicalDevices/ResourcesforYou/Industry/default.htm>.

Sincerely yours,



for Anthony D. Watson, B.S., M.S., M.B.A.
Director
Division of Anesthesiology, General Hospital,
Infection Control and Dental Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

Indications for Use

510(k) Number (if known): **K113520**

Device Name: **SYSTEM 1E Liquid Chemical Sterilant Processing System**

Indications For Use:

The SYSTEM 1E Liquid Chemical Sterilant Processing System is intended for liquid chemical sterilization of cleaned, immersible, and reusable critical and semi-critical heat-sensitive medical devices in healthcare facilities.

The SYSTEM 1E Processor dilutes the S40 Sterilant Concentrate to its use dilution (>1820 mg/L peracetic acid), liquid chemically sterilizes the load during a controlled 6-minute exposure at 46 to 55°C, and rinses the load with extensively treated* potable water. The processed load should be used immediately.

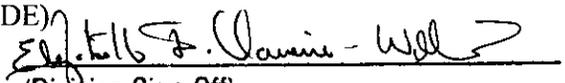
The SYSTEM 1E Processor uses only S40 Sterilant Concentrate to liquid chemically sterilize medical devices.

- * The extensive treatment of EPA potable water consists of:
1. Pre-filtration through two pre-filters:
 - Pre-filter 1 is a gross depth filter that removes approximately 5 micron or larger particles/contaminants.
 - Pre-filter 2 is a surface filter that removes particles/contaminants > 0.1 micron.
 2. UV Irradiation:
 - During transit through the UV water treatment chamber, a UV dose sufficient to achieve a ≥ 6 -log reduction of MS2 virus is delivered to the water.
 3. 0.1 micron filtration:
 - The water prepared by pre-filtration and UV irradiation is filtered through redundant, 0.1-micron (absolute rated) membranes to remove bacteria, fungi and protozoa > 0.1 micron.

Prescription Use _____ AND/OR Over-The-Counter Use X
(Part 21 CFR 801 Subpart D) (21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)


(Division Sign-Off)
Division of Anesthesiology, General Hospital
Infection Control, Dental Devices