



QUALITY FOR LIFE

K070318

MAY 11 2007

510(k) SUMMARY  
of  
SAFETY and EFFECTIVENESS

**A. General Information**

1. *Submitter's Name:* OTTO BOCK HealthCare LP
2. *Address:* Two Carlson Parkway N., Suite 100  
Minneapolis, MN 55447-4467
3. *Telephone:* 763-253-5607
4. *Contact Person:* Shannon Tufvander
5. *Date Prepared:* January 18, 2007
6. *Registration Number:* 2182293

**B. Device**

1. *Name:* Skippi Powered Wheelchair for Children
2. *Trade Name:* Skippi Powered Wheelchair for Children
3. *Common Name:* Powered wheelchair
4. *Classification Name:* Powered wheelchair
5. *Product Code:* ITI
6. *Class:* II
7. *Regulation Number:* 890.3860

**Customer Service, Technical Support, and General:** 800.328.4058 • Fax 800.962.2549

**Fabrication Orders:** 877.FAB.OTTO (877.322.6886) • Fax 800.810.7994

**Customer Satisfaction Hotline:** 877.627.6583

[www.ottobockus.com](http://www.ottobockus.com)



QUALITY FOR LIFE

### C. Identification of Legally Marketed Devices

1. *Name:* A-200
2. *K Number:* K052706
3. *Date Cleared:* December 27, 2005

### D. Description of the Device

The Skippi is propelled by two DC Motors, Micro Motor AG PMGST 76/40 L/R 6 Km/H at 24 Volts. The controller is the enAble 40 manufactured by Curtis Instruments, AG.

The closed "H" frame is rigid (not foldable), all four wheels are suspended independently. The Skippi has 12.5 inch knobby rear tires.

Options available to the end user are the following:

- Seat Tilt Mechanism, Electric or with Gas Compression Spring
- Electric Backrest Angle Adjustment
- Seating Shell Adapter
- Tray
- Control for Attendant
- Bumper Bar
- Single Panel Footrest
- Colors (Yellow, Blue, or Red)

### E. Intended Use Statement

The Skippi is a rear wheel drive powered wheelchair with caster front wheels for indoor use (primarily) but can be used outdoors as well. These wheelchairs provide mobility to physically challenged persons. The wheelchair can be moved by the user operating the Curtis Instruments enAble 40 that is connected to the Micro Motor. The wheelchair is steered by different rotation of the rear wheels.

### F. Technological Characteristics Summary



**QUALITY FOR LIFE**

The Skippi Wheelchair is substantially equivalent to the OTTO BOCK A-200 Wheelchair, cleared on December 27, 2005 as K052706.

Each wheelchair is a powered wheelchair for the active user, with a rigid frame and similar characteristics. The Skippi is for children or small people up to 110 pounds, whereas the A-200 is for anyone up to 200 pounds.

The Skippi was tested by TÜV Product Service to the following standards:

- EN 12184
- ISO 7176 -- Series
- ANSI/RESNA WA Vol. 2 Section 21 Amendments 1998 for EMC

with the conclusion that “the test sample fulfills the requirements.”



Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

OTTO BOCK HealthCare, LP  
% Mr. William Jackson  
Regulatory Consultant  
2247 Jennifer Lane  
St. Paul, Minnesota 55109-2851

MAY 11 2007

Re: K070318

Trade/Device Name: Skippi Powered Wheelchair for Children  
Regulation Number: 21 CFR 890.3860  
Regulation Name: Powered wheelchair  
Regulatory Class: II  
Product Code: ITI  
Dated: April 18, 2007  
Received: April 27, 2007

Dear Mr. Jackson:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

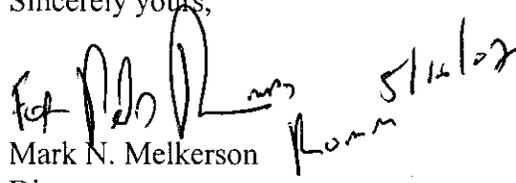
Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

Page 2 – Mr. William Jackson

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (240) 276-0120. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (240) 276-3150 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

A handwritten signature in black ink, appearing to read "Mark N. Melkerson", with a date "5/12/02" written to the right of the signature.

Mark N. Melkerson  
Director  
Division of General, Restorative  
and Neurological Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

**Indications for Use**

**510(k) Number** (if known): *To be determined*

**Device Name:** Skippi Powered Wheelchair for Children

**Indications for Use:**

- Provide mobility to persons physically challenged and limited to sitting positions.

Prescription Use \_\_\_\_\_  
(Part 21 CFR 801 Subpart D)

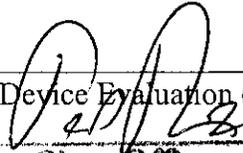
AND/OR

Over-The-Counter Use  X   
(21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE –  
CONTINUE ON ANOTHER PAGE IF NEEDED)

---

Concurrence of CDRH Office of Device Evaluation (ODE)

  
\_\_\_\_\_  
(Division Sign-Off)

**Division of General, Restorative,  
and Neurological Devices**

**510(k) Number** 16070318