



JUN 16 2008

Suite 316
1800 Sandy Plains Industrial Parkway
Marietta, GA 30066-6363
Phone: 770.425.4470 Fax: 770.425.4740
www.chart-ind.com

510(K) Summary

K080023

Submitter: Caire, Inc
1800 Sandy Plains Industrial Parkway
Marietta, GA 30066-6363

Contact Person: Zulfiqar Ali Choudary
BioMedical Quality & Regulatory Affairs Manager

Phone Number: 770-257-1276
Facsimile Number: 770-425-4740
E-Mail Address: ali.choudary@chart-ind.com

Date Prepared: November 12, 2007

Device Name: Spirit 300
Classification Name: Portable Liquid Oxygen Unit
Medical Device Class: II
Classification Panel: Anesthesiology
Product Code: BYJ
CFR Section: 865.5655
Regulation Description: Portable Liquid Oxygen Unit

Predicate Device: Spirit 300 [FDA 510(K) K013251]

Device Description:

The spirit 300 is a small, lightweight, liquid oxygen portable unit. The unit consists of a vacuum insulated cryogenic dewar, a vaporizer coil, an economizer regulator, two pressure relief valves, a manifold, an electronic printed circuit board conserving device (digital) and a protective case. The dewar has a capacity and capability to store 0.3L of liquid oxygen. The economizer regulator either draws gas from the dewar head space, if the internal pressure is above 20 PSI. If the internal pressure is below the 20 PSI, the economizer regulator will draw liquid oxygen, from the bottom of the dewar, which must then have to pass through a long vaporizing coil. The economizer regulator assures a

longer duration time for the unit, because it allows the gas in the dewar headspace to be used for oxygen delivery. The gas outlet of the economizer regulator connects with the manifold, dependent on the flow selection setting, the manifold either delivers the gas to the conserving device PCB or diverts a 2L PM continuous flow directly to the oxygen outlet. The continuous flow setting is labeled CF on the flow selector. The conserving device has equivalent settings for 1, 1.5, 2, 3, 4 and 5 LPM prescriptions. At the various conserving device settings, the outlet gas is delivered in quick pulse dosage just at the onset of patient inhalation. The amount of gas delivered (with each breath) is approximately 15 ml/LPM setting.

Indications for use:

The Caire Spirit 300 will provide a source of supplemental oxygen for ambulatory healthcare patients by vaporizing medical grade liquid oxygen and then dispensing it to the patient via an integral digital electronic conserving device (DOCD). The spirit 300 is neither a life sustaining nor life supporting device.

Statement of Technical Characteristics Comparison:

Caire Inc confirms that the modifications described here within have not altered the fundamental scientific technology of the predicate device (Spirit 300) cleared under 510(K) number K013251.



Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

JUN 16 2008

Mr. Zulfiqar Ali Choudary
Quality and Regulatory Affairs Manager
Caire, Incorporated
1800 Sandy Plains Industrial Parkway, Suite 316
Marietta, Georgia 30066

Re: K080023
Trade/Device Name: Spirit 300
Regulation Number: 21 CFR 868.5655
Regulation Name: Portable Liquid Oxygen Unit
Regulatory Class: II
Product Code: BYJ
Dated: May 21, 2008
Received: May 27, 2008

Dear Mr. Choudary:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.


If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (240) 276-0120. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

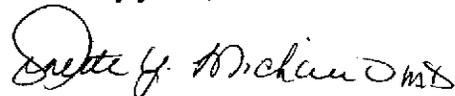

Chiu Lin, Ph.D.

Director

Division of Anesthesiology, General Hospital,

Infection Control and Dental Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

Indications for Use

510(k) Number (if known):

Device Name: Spirit 300

Indications for Use:

The Caire Spirit 300 will provide a source of supplemental oxygen for ambulatory healthcare patients by vaporizing medical grade liquid oxygen and then dispensing it to the patient via an integral electronic conserving device (digital oxygen conserving device-DOCD). The Spirit 300 is neither a life sustaining nor life supporting device.

Prescription Use X
(Part 21 CFR 801 Subpart D)


and/or

Over-the-counter Use _____
(21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE, CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Page 1 of 1

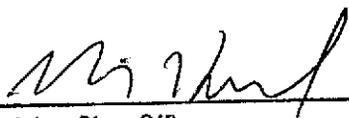

(Division Sign-Off)
Division of Anesthesiology, General Hospital
Infection Control, Dental Devices

510(k) Number: K080023